

# Tryton Side Branch Stent

**TRYTON**<sup>®</sup>  
Medical  

---

BUILT FOR BIFURCATION




New SHORT Stent Length







# Dedicated Bifurcation Stent

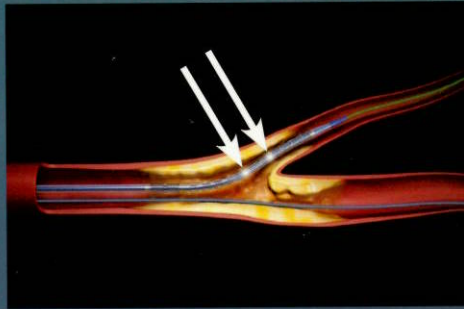

Position Side Branch Stent

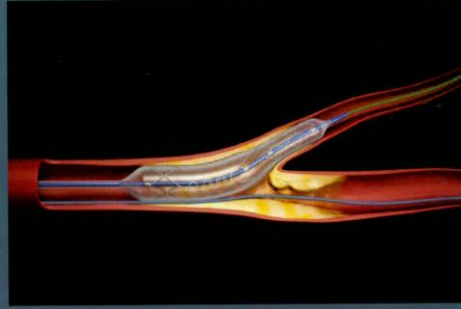


Deploy Side Branch Stent

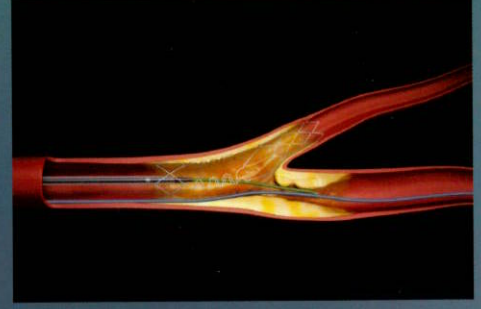

Advance Wire in Main Branch



Position Main Vessel DES

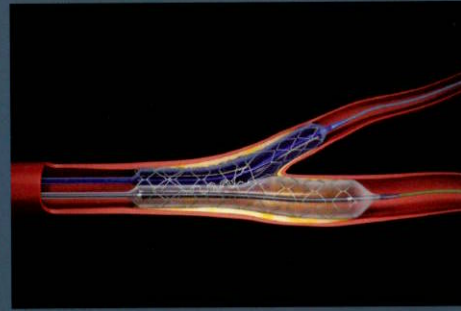

Perform Simultaneous Kissing Balloon Post Dilatation



The Bifurcation is Protected and Secured

Model	Reference Number	Diameter SB - MB (mm)	Length (mm)	Minimum Guiding Catheter Diameter	Maximum Post-expansion Diameter SB - MB (mm)	Nominal Pressure (atm)	RBP (atm)
<b>STANDARD</b>							
Straight	T5-2525-191	2.5 - 2.5	19	5F	3.0 - 4.0	8	14
Step	T5-2530-191	2.5 - 3.0	19	5F	3.0 - 4.0	10	14
Step	T5-2535-191	2.5 - 3.5	19	5F	3.0 - 4.0	10	14
Step	T5-3035-181	3.0 - 3.5	18	6F	4.0 - 4.5	10	14
Step	T5-3540-181	3.5 - 4.0	18	6F	4.0 - 4.5	8	14
<b>SHORT</b>							
Step	T5-3035-151	3.0 - 3.5	15	6F	4.0 - 4.5	8	14
Step	T5-3540-151	3.5 - 4.0	15	6F	4.0 - 4.5	10	14


Investigational in the US and not approved in Japan. The Tryton Side Branch Stent System is for sale in Europe.

**TRYTON**  
Medical  
BUILT FOR BIFURCATION

Corporate Headquarters  
1000 Park Forty Plaza, Suite 325  
Durham, NC 27713 USA  
Phone: +1.919.226.1490  
Fax: +1.919.226.1497

Tryton Medical B.V.  
Centaurusweg 123  
5015TC Tilburg  
The Netherlands

Distributed by

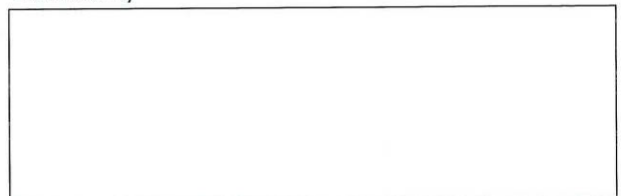

For more information please contact us at [www.trytonmedical.com](http://www.trytonmedical.com)

©2013 Tryton Medical, Inc. All rights reserved. D2016G\_V06